


RESUSCITATION / VENTILATION

Ventlab

AIRFLOW™ BAG - SS3200MB

LATEX FREE | SINGLE PATIENT USE | DISPOSABLE | CE MARK | DEHP

- Single use resuscitation bag
- Infant face mask
- Non-sterile and latex free
- Volume - 250ml
- Stroke Volume - 140-170ml
- Weight Range - ≤10kg (22lbs)

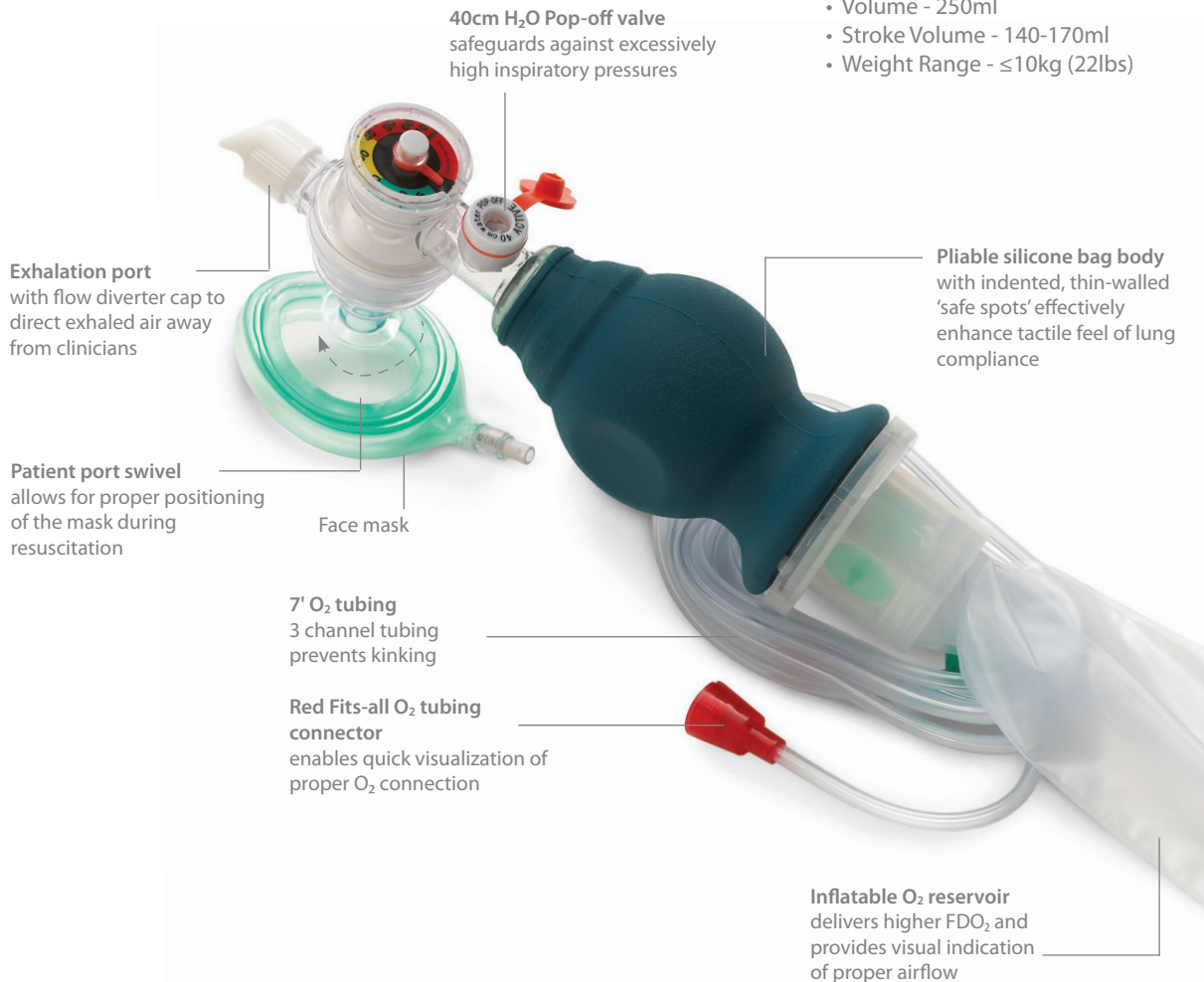

INSTRUCTIONS FOR USE

PREPARATION FOR USE

1. Test the resuscitator for proper functioning: With patient port completely occluded, squeeze bag body to assure resistance is present and positive needle movement occurs.
2. Before using the mask, inspect for adequate inflation.
3. The manual resuscitator/ventilator may be used with a 19 mm ID PEEP accessory. Attach PEEP accessory to the exhalation port. Be sure that the accessory fits properly and does not interfere with compression of the resuscitator.
4. Actual PEEP may vary with patient lung compliance and resistance. Verify PEEP with a certified manometer.
5. For correct performance on the manual resuscitator/ventilator with oxygen reservoir, unfold the reservoir bag and assure that airflow is not restricted.

CAUTIONS

On models with option Pop-Off feature: To override the Pop-Off feature, insert tethered cap into Pop-Off opening.

DIRECTIONS FOR USE

1. Place the patient in a supine position. Establish and maintain an open airway.
2. Grasp the bag body with one hand. If a more sensitive indication of lung compliance is desired, grasp bag body at the two "safe spot" areas.
3. Hold the mask between the index finger and thumb of the other hand. Place mask over face firmly to form a tight seal around the patient's nose and mouth.
4. Ventilate the patient by compressing the bag body for inhalation and releasing the bag body for the patient's passive exhalation and bag body re-expansion. Continue this cycle as directed by medical authority.
5. Monitor peak airway pressure by observing the built-in gauge.
6. To remove vomitus: Disconnect resuscitator from patient. Tap the patient valve several times while squeezing the bag body. Re-test the resuscitator for proper functioning.
7. If patient is intubated, remove mask from patient port. Connect patient port directly to the endotracheal tube adapter. Continue ventilation.
8. To use supplemental oxygen: Connect oxygen supply tubing to O₂ source at appropriate flow rate. FDO₂ values may be affected if flow is not sufficient. Oxygen flow ≥15 lpm may be necessary. **Do not let flow rate exceed 30 lpm due to possible increase in exhalation resistance.**
9. Replace resuscitation bag when it is visibly soiled or per hospital policy — whichever comes first

WARNING: Always read instructions for use on product packaging prior to use on of these devices.

