


OXYGEN DELIVERY

Ventlab

OXYGEN MASKS

LATEX FREE | SINGLE PATIENT USE | DISPOSABLE | CE MARK

The Medium Concentration Oxygen mask delivers 6-10 liters per minute with an oxygen concentration up to 50%. All masks include 7' 3-Channel oxygen tubing. Additional clear or semi-transparent green supply tubing is available in 6 lengths.

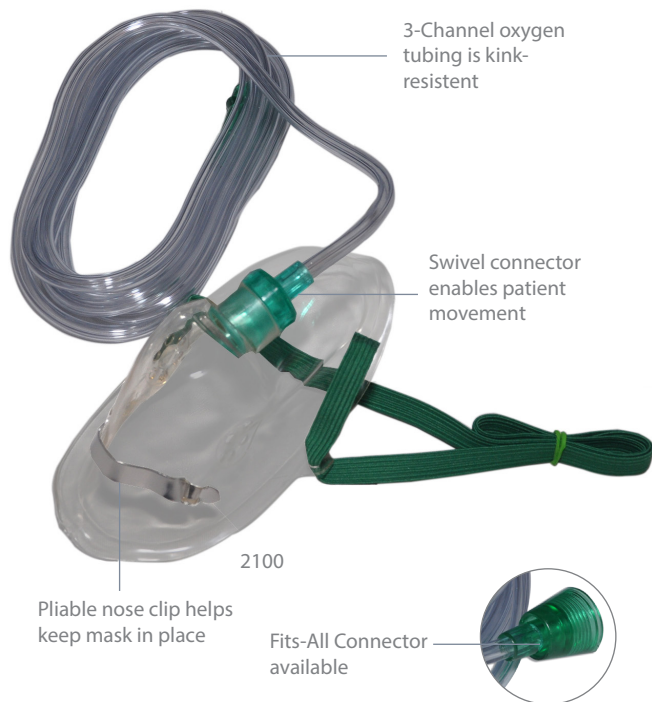

MEDIUM CONCENTRATION

- Flow rate of 6-10 liters per minute
- Delivers an oxygen concentration of 40% to 50% based on breathing rate and tidal volume
- Soft, pliable mask fits comfortably on patient
- Includes elastic head strap for best fit
- Transparent PVC mask enables clear visualization of patient
- O₂ Inlet designed to direct oxygen into the nostrils rather than up toward eyes
- All sizes include 7' 3-Channel oxygen tubing
- Pliable nose clip helps keep mask in place
- Swivel connector enables patient rotation and movement without discomfort

OXYGEN SUPPLY TUBING

- Available in clear or semi-transparent green
- Six lengths available
- Three-channel tubing design prevents kinking

CAT#	DESCRIPTION	SIZE	PK
2100	Medium Concentration	Adult	50
2100F	Medium Concentration + Fits-All Connector	Adult	50
2200	Medium Concentration	Pediatric	50
2200F	Medium Concentration + Fits-All Connector	Pediatric	50
2300	Medium Concentration	Infant	50
2300F	Medium Concentration + Fits-All Connector	Infant	50

CAT#	DESCRIPTION	TUBING LENGTH	PK
3004	Clear Supply Tubing	4'	50
3007	Clear Supply Tubing	7'	50
3007G	Green Supply Tubing	7'	50
3007F	Clear Tubing + Fits-All Connector	7'	50
3015	Clear Supply Tubing	15'	50
3015G	Green Supply Tubing	15'	50
3025	Clear Supply Tubing	25'	50
3025G	Green Supply Tubing	25'	50
3040	Clear Supply Tubing	40'	50
3040G	Green Supply Tubing	40'	50
3050	Clear Supply Tubing	50'	50
3050G	Green Supply Tubing	50'	50



Ventlab OXYGEN MASKS

LATEX FREE | SINGLE PATIENT USE | DISPOSABLE | CE MARK

Non-rebreathing masks deliver a high concentration of oxygen with a partial non-rebreather option available. Non-rebreathers feature a swivel connector for optimum patient comfort and uninterrupted oxygen flow.



HIGH CONCENTRATION OXYGEN MASKS

- Delivers high concentration of oxygen
- Reservoir bag assures oxygen supply to meet variable breathing patterns and tidal volumes
- Soft, pliable mask fits comfortably on patient
- Includes elastic head strap for best fit
- Transparent PVC mask enables clear visualization of patient
- Swivel connector design ensures uninterrupted flow of oxygen from the bag to the patient
- All sizes include 7' 3-channel oxygen tubing
- Pliable nose clip helps keep mask in place
- Available as non-rebreather or partial non-rebreather
- Designed with two valves on the non-rebreather; one on the partial non-rebreather

PARTIAL NON-REBREATHER ONE VALVE

CAT#	DESCRIPTION	SIZE	PK
2101	Partial Non-Rebreathing	Adult	50
2101F	Partial Non-Rebreathing + Fits-All Connector	Adult	50
2201	Partial Non-Rebreathing	Pediatric	50
2201F	Partial Non-Rebreathing + Fits-All Connector	Pediatric	50
2301	Partial Non-Rebreathing	Infant	50
2301F	Partial Non-Rebreathing + Fits-All Connector	Infant	50

NON-REBREATHER TWO VALVES

CAT#	DESCRIPTION	SIZE	PK
2102	Non-Rebreathing	Adult	50
2102F	Non-Rebreathing + Fits-All Connector	Adult	50
2202	Pediatric	Pediatric	50
2202F	Non-Rebreathing + Fits-All Connector	Pediatric	50
2302	Infant	Infant	50
2302F	Non-Rebreathing + Fits-All Connector	Infant	50

