



CONTRACT NUMBER HPG-2543

August 1, 2015 - September 30, 2022

Infu-Surg™

Infu-Surg™ Clarity

Available with
economical **mesh** panel
or **clear** panel



Ethox

PRESSURE INFUSION BAGS


- Pressure relief valve prevents over inflation (330mmHg pressure relief)
- Large, oval shaped bulb allows for quick and easy inflation of bladder
- Dedicated fluid bag hook enables efficient loading and unloading of the system without having to remove the infusion bag from the IV pole
- Suitable for use with external inflation sources
- Color-coded gauge makes for accurate pressure monitoring (0-300mmHg)
- Incredibly reliable – 100% tested



Clear solution bag panel enables patient label and barcode scanning without having to remove bag.



Three-way stopcock ensures precise control of pressure while the large oval bulb allows for quick bladder inflation, all with ONE hand.

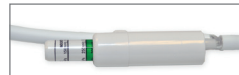

360° view of the color-coded gauge facilitates accurate pressure monitoring while the safety pressure relief valve prevents over-inflation.

Infu-Surg™

500ml 4005H
1000ml 4010H
3000ml 4030



Infu-Surg™ Clarity

500ml 4050
1000ml 4100



Product details available
at www.Sun-Med.com

For more information,
call **1-800-433-2797**