


COMBIGUARD II® + COMBIPORT®

LATEX FREE | STERILE | SINGLE PATIENT USE | DISPOSABLE

The Combiguard II® and Combiport® are easy to use, highly effective wound irrigation tools. The Combiguard II® Irrigation Splash Guard covers the wound and helps prevent splash-back. The Combiport® attaches to a standard saline bag to deliver a high volume of fluid rapidly. Available sterile and non-sterile.

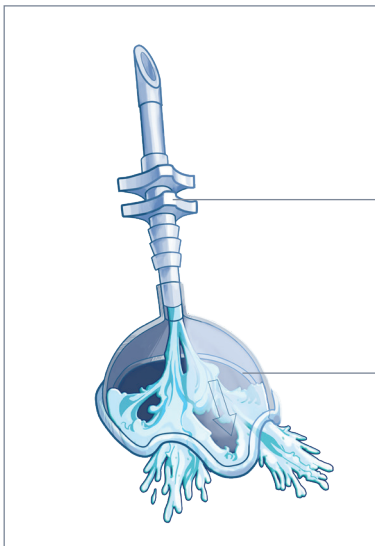

Working together, a saline bag uses a Combiport attached to a Combiguard II to effectively cleanse the wound without splashback on the caregiver

The Combiguard II attaches to a syringe (or Combiport) to assist in wound cleansing

Fluid flows through angled port, directed away from cleansed area and toward the outlet ports



COMBIPORT®
Wound Irrigation Device

COMBIGUARD II®
Irrigation Splash Guard

COMBIPORT® WOUND IRRIGATION DEVICE

- Attaches to a standard saline bag to administer a high volume of fluid in a short period of time
- Irrigation helps remove foreign bodies, reduce pathogen numbers and moisten wound tissue
- Patented spike attaches to standard 500 mL or 1000 mL saline bag
- For use only with Combiguard Irrigation Splash Guard (not sold separately)
- Removable cap for easy stop/start during cleansing
- Easy to assemble and use

COMBIGUARD II® IRRIGATION SPLASH GUARD

- Egg-shaped cup covers wound site to prevent irrigation splash back, helping to protect clinician
- Dual outlet ports enable contaminated fluids to safely escape
- Fluid delivery inlet delivers a stream equivalent to an 18 g needle
- Angled delivery inlet directs fluid away from cleansed area toward outlet ports
- Rounded edges provide a smooth surface against the skin
- Fits on 30-40 cc luer-tipped syringe
- Can be combined with Combiport Wound Irrigation Device

STERILE

ITEM	DESCRIPTION	PK
10222	Combiguard II® Irrigation Splash Guard	50
10100	Combiguard II® + Combiport®	50

NON-STERILE

ITEM	DESCRIPTION	PK
10222NS	Combiguard II® Irrigation Splash Guard Non-Sterile	500